

Distributed by:



Carl Stuart Limited

ADVANCED APPLIED TECHNOLOGIES

Contact Us:

Irl Ph: 01 4523432

UK Ph: 08452 30 40 30

Web: www.carlstuart.com

Email: info@carlstuart.com

**Operating Instructions for the
Disk Mill Model DM 200**

Retsch®

Information on these Operating Instructions

The present operating instructions for the laboratory model DM 200 disk mill provides you with all of the necessary information on the areas specified in the table of contents.

It is an instruction for the target group(s) defined for each of these areas in the safe and intended use of the **DM 200**. The prerequisite of each of the target group(s) being able to use it safely and as intended is having a knowledge of the relevant sections.

The present technical documentation is a reference work and instructions for learning and each of the sections is closed in and of itself.

These operating instructions are not instructions for repair. If repairs should be necessary, please consult your vendor or get in contact with Retsch GmbH directly.

Information on these Operating Instructions	2
Safety	4
Safety Information.....	4
Warning Signs	4
Repairs	5
Confirmation	5
Technical Data.....	6
Use as intended.....	6
Functional Description	6
Protective Equipment.....	7
Drive	7
Motor Output	7
Rated Voltages.....	7
Rpm's.....	7
Protective Types	7
Emissions	8
Machine Dimensions.....	8
Floor Space Required	8
Transportation and Installation	9
Packaging.....	9
Transportation.....	9
Intermediate Storage	9
Scope of Delivery.....	9
The Parameters for the Place of Installation	10
Installation	10
Electrical Connection	10
Operation	11
The Operating Elements and Operation	11
Adjusting the Gap Width.....	12
Starting the Machine	12
Feeding the Material to be Crushed.....	12
Changing the Direction of Rotation and Replacing the Milling Disc	13
General Information.....	14
Accessories	14
Cleaning.....	14
Service.....	14
Necessary Safety Tests	14
Copyright	14
Modifications	14
Safety Regulations (Table).....	15
The Terms of the Warranty	16

Safety

The DM 200 laboratory disk mill is a highly modern high-performance product from Retsch GmbH and is state-of-the-art. If the machine is used as intended and if the user has a knowledge of the technical documentation here, it is completely safe when operated.

Safety Information

You as the operator have to ensure that the persons entrusted with working at the **DM 200**:

- have taken notice of and understood all of the regulations in the area of safety,
- know all of the instructions for use and regulations of the target group relevant for them before commencing work,
- have access to the technical documentation of this machine at any time and without any problems.

Before work is commenced on the **DM 200**, ensure that new personnel has been made familiar with its safe and intended use either by having a competent person introduce it to them orally and / or based upon the present technical documentation.

Improper operation may lead to personal injury, property damage or injury. You are responsible for your own safety and that of your employees.

Ensure that no unauthorised persons have access to the DM 200 laboratory disk mill.

For your own protection, have your employees confirm that they have been instructed in the operation of the **DM 200**. You will find the draft of the appropriate form after the section on safety.



We will not be liable for any claims to compensation whatsoever for property damage or personal injury that were caused by not observing the following safety information.

Warning Signs

we give warnings with the following signs:



personal injury



property damage

Repairs

These operating instructions are not instructions for repair. For your own safety, repairs may only be done by Retsch GmbH or an authorised agency or service technicians.

In this case, please notify:

The Retsch agency in your country
your vendor
Retsch GmbH directly

your service address:

--

Confirmation

I have taken notice of the foreword and the section on safety.
 <hr/>
the signature of the operators
 <hr/>
the signature of the service technician

Technical Data

machine model designation: **DM 200**

Use as intended

The **DM 200** laboratory disk mill is suited for finely milling soft to hard materials (up to 8 Mohs) in sets or continuously from the fields of mining and metallurgy, rocks and earth, the glass industry or soil research.

Our applications laboratory would be glad to give you further information upon request.

The **DM 200** is designed for quantities from 20 to approximately 150 Kg/h depending upon the setting of the discharge gap and the piled density and crushing behaviour of the sample. The feed size may not exceed an edge length of a maximum of 20 mm.



The DM 200 is not designed to be explosion protected and therefore it is not suited to milling explosive self-igniting or fire-promoting substances.



If oxidised materials are milled such as metals, organic materials, wood, plastics, etc., there is the danger of self-ignition (dust explosion) if the proportion of fines exceeds a certain percentage. It is necessary to observe the appropriate safety guidelines when materials such as these are crushed.



Do not modify the machine and only use the spare parts and accessories approved by Retsch.

Otherwise the conformity with European directives declared by Retsch shall be rendered invalid.

Furthermore, this leads to a loss of all warranty claims.




Functional Description

The material sample is crushed by two milling discs coarsely meshed on the inside working against one another in a milling space sealed against dust from the outside. One of these is driven by a powerful and slow moving gear motor.

The material to be milled is filled into the centre of the stationary milling disc through a closable funnel and it escapes from the gap between the two disks after being crushed by compressing and shearing force. The width of the gap can be adjusted and controlled from outside and determines the mean granulation of the material to be milled.

Protective Equipment

The following safety equipment reliably prevents persons from unintentionally putting their hands in the machine when it is running.

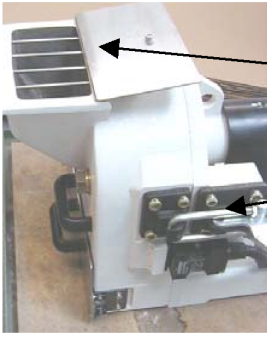


- funnel tube designed in conformity with EN 294.
This prevents persons from putting their hands into the milling space when the flap is opened

- safety limit switch.
This prevents the machine from starting up when the milling space is open.
It also stops the motor within a safe period of time when the milling space is opened during operation



- safety limit switch.
This prevents the machine from starting up when the collecting container is drawn out.
It also stops the motor within a safe period of time when drawing out the collecting container.

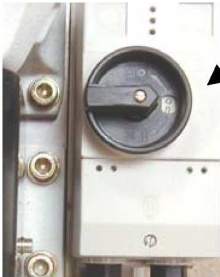

- protective motor switch, ON/OFF switch.
This switches off the machine automatically if there is a defect or if the motor is overloaded.
position [I] the machine starts.
position [0] the machine stops.
This prevents the machine from starting up independently if there is an electrical defect or power failure.

Drive

standard three-phase gear motor

Motor Output

1,500 Watt

Rated Voltages

3~ 230 V	50 / 60Hz	approx. 5.7 A
3/N~ 400 V	50 / 60Hz	approx. 3.3 A
3~115 V	50 / 60Hz	approx. 11.4 A
3~ 200 V	50 / 60Hz	approx. 6.5 A

Rpm's

motor Rpm's 2,850 min⁻¹

Protective Types

gear motor	IP 55
safety limit switch for the door	IP 67
safety limit switch for the collecting container	IP 67
protective motor switch, ON-OFF switch	IP 55

Emissions

Characteristic Noise Values :

Noise measurement in conformity with DIN 45635-031-01-KL3

The properties of the material to be crushed also has an effect upon the characteristic noise values.

sound level $L_{WA} = 81 \text{ dB(A)}$

emission read at the workpiece $L_{pAeq} = 69,4 \text{ dB(A)}$

Machine Dimensions

height: 400 mm, 500 mm with the vertical flap on the funnel.

width: 430 mm, 900 mm with door opened at 180°.

depth : 890 mm, 1,000 mm with door opened at 90°

weight: approximately 140 kg

Floor Space Required

430 mm x 1000 mm; a safety distance is not necessary.

Transportation and Installation

Packaging

Packaging is in accordance with the transportation path and satisfies the generally recognised packaging guidelines.



Please store the packaging for the duration of the warranty period because your warranty claim would be endangered if there were a complaint and it were returned in insufficient packaging.



Transportation

In order to be able to transport the **DM 200** in a crane of a minimum load-bearing capacity of 250 kp, it should be hoisted with the hoisting equipment as shown in the adjacent figure.



The DM 200 may not be shocked, shaken or thrown during transportation. Otherwise, the electronic and mechanical component parts could be damaged.



Temperature Fluctuations

If there are strong temperature fluctuations (for instance, when being transported by aeroplane), the DM 200 should be protected from condensation water. Otherwise, there could be damage to the electronic component parts.

Intermediate Storage

Please also ensure that the DM 200 is stored dry if it is placed in intermediate storage.

Scope of Delivery

- the DM 200 disk mill
- the operating instructions

Please check to see whether the delivery is complete including the accessories specifically ordered.

Check to see whether the **DM 200** is in perfect working order (refer to the section on operation for this)



If the delivery is incomplete and / or if there was damage during transportation, you have to notify the transporter and Retsch GmbH without delay (within 24h). Under certain circumstances, later complaints may not be recognised.

The Parameters for the Place of Installation

ambient temperature :

5°C to 40°C



If the temperature is above or below the ambient temperature, the electronic and mechanical component parts may be damaged changing the performance data in an unpredictable fashion.

humidity :

The maximum relative humidity is 80% at temperatures up to and including 31°C, decreasing linearly to 50% relative humidity at 40°C



If the humidity is higher, the electronic and mechanical component parts may be damaged changing the performance data in an unpredictable fashion.

Installation Altitude:

a maximum of 2,000 m above sea level

Installation

Only place the DM 200 on a solid and sturdy foundation.

Electrical Connection

The mains fusing should be designed as follows:

at 3~ 230 and 3/N~ 400V = 16A

at 3~ 115 and 3~ 200V = 20A

- you may find the voltage and frequency of the **DM 200** on the nameplate.
- please ensure that the values agree with the existing power mains.
- only connect the DM 200 to the power mains with the aid of the cable connector supplied.

You may have the connection made directly with a change-over switch to avoid having to call in an electrician. This has the benefit that the milling disc teeth may be used by both sides. This is described in the section "**Changing the Direction of Rotation and Replacing the Milling Disc**".



Electrical and mechanical component parts may be damaged if the values on the nameplate are not complied with.





Modifications or rebuilding it to a supply voltage other than stated on the machine nameplate should only be done by an electrician.

Danger of Current Surge!

Operation

The Operating Elements and Operation

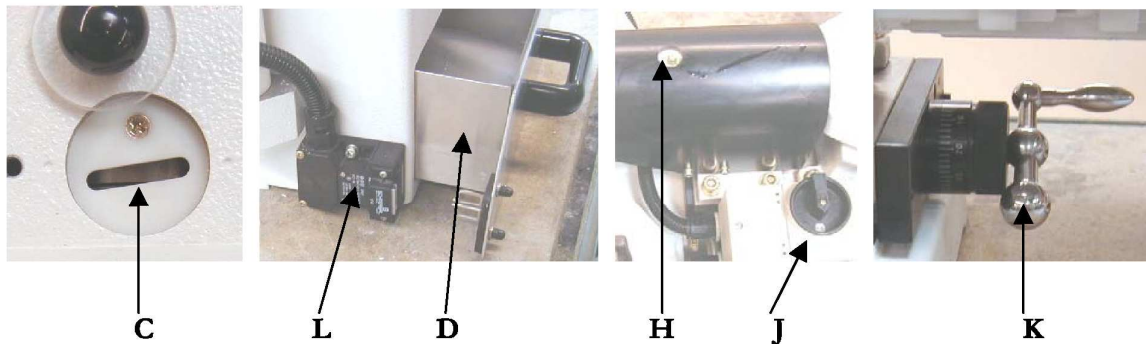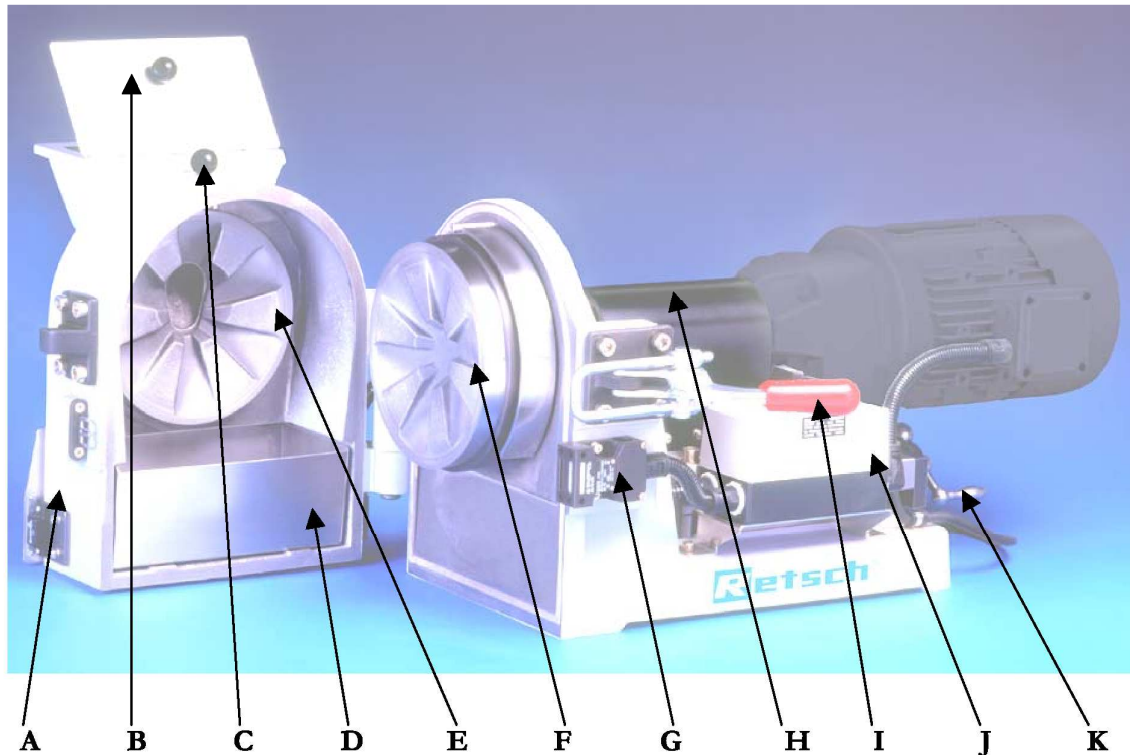


item	element	function
A	door	seals the milling space dustproof
B	funnel tube	receives the material to be milled
C	window milling gap	enables you to check the milling gap by means of blade gauge
D	collecting container	receives the crushed material to be milled
E	milling disc doors	crushes the material to be milled fed into it together with F
F	milling disc housing	crushes the material to be milled fed into it together with E
G	safety switch for the door	prevents the machine from starting up when the door is open
H	lubrication point	enables you to relubricate the coupling
I	door lock	enables you to open and safely close the doors
J	protective motor switch	disconnect the motor from the mains if there is an overload position [I] the machine is switched on position [0] the machine is stopped
K	gap adjustment crank	enables you to adjust the gap together with C graduation line = 0.01 mm
L	safety switch for the collecting container	prevents the machine from starting up when the collecting container is drawn out

Adjusting the Gap Width

You may adjust the gap width between the milling disc in the housing and in the door just by adjusting the graduated rotating crank mounted below the motor. 1 graduation line = 0.01 mm.


You may adjust the gap from 0.1 to approximately 5 mm.

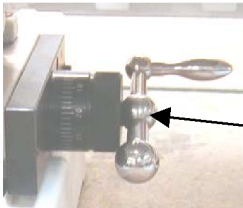

- **Rotating Crank**

Caution! Do not go below the minimum gap width of 0.1 mm.

You may control check the gap width setting through a window behind the funnel tube by means of a sensing gauge.

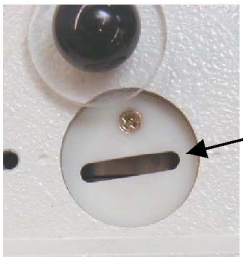

- **Window**

You may also adjust the gap width when the machine is running.

Material that is difficult to mill should be milled at least twice. Select a larger gap width as pre-crushing for the first pass and then you achieve the desired final fineness with the small gap in the second milling pass.



Caution !


The smallest gap width is 0.1 mm.

The milling discs may not touch one another because otherwise this could damage the milling discs.

Starting the Machine

It is only possible to start the DM 200 when the door is closed and the collecting container has been inserted.

The ON-OFF switch is on the right-hand side of the **DM 200**.

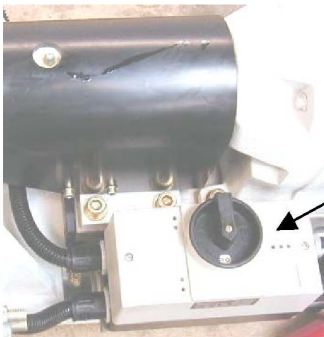

- **ON-OFF switch**

position [I] **DM 200** starts

position [0] **DM 200** stops

The ON-OFF switch also acts as a protective motor switch.


If the motor is overloaded or if there is another electrical fault, this switch disconnects the **DM 200** from the mains.

The **DM 200** can only start up when the **DM 200** is manually started with the ON-OFF switch.

Feeding the Material to be Crushed

After the gap has been set and the **DM 200** has been started, you commence material feed. Do not exceed a maximum edge length 20 mm.

Open the lid of the funnel tube for this purpose.

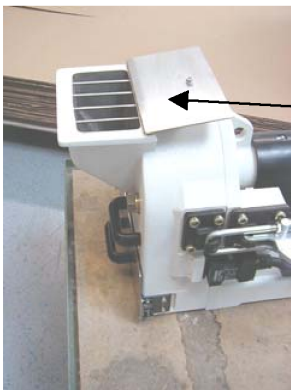

- **the lid of the funnel tube**

When feeding the material to be crushed, please ensure that the funnel tube lid is closed because otherwise the material to be crushed could fall back.

Furthermore, the quantity of fed material to be crushed depends upon its millability. Therefore, note the decreasing intensity of the milling sound to ascertain the optimum feed quantity.



Caution!


Only feed the material to be crushed after having started the **DM 200**.

If the **DM 200** starts up filled with material to be crushed, this may lead to damage in the mechanical component parts.

Changing the Direction of Rotation and Replacing the Milling Disc

The milling discs are subject to natural wear and tear after a longer period of time of use.

Before they have to be replaced with new ones, you may change the direction of the motor's rotation to be able to use the opposite side of the gear teeth. This also extended the milling disc service life.

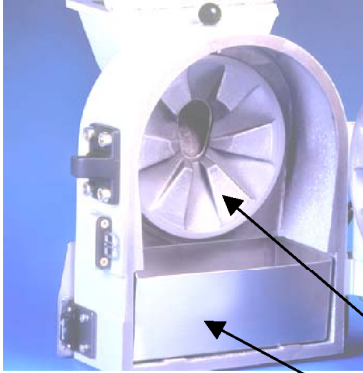


Changing the Direction of Rotation :

The direction of the motor's rotation may be changed by changing the phases on the network. Since this requires interfering with the electrical power leads to the DM 200, only an electrician may change the direction of the motor's rotation.

When it is connected up to the network for the first time, you may have a change-over switch mounted to avoid having to call in an electrician if the direction of the motor's rotation is changed again.


Replace the Milling Disc on the Door :

- pull the mains plug
- remove the collecting container
- keep the door closed
- loosen the hexagonal screws slightly with a SW 30 open-jawed spanner
- open the door
- hold the milling disc securely and unscrew both hexagonal screws
- remove the milling disc
- clean the receptacle of milling disc in the door
- insert a new milling disc
- insert the hexagonal screws again
-

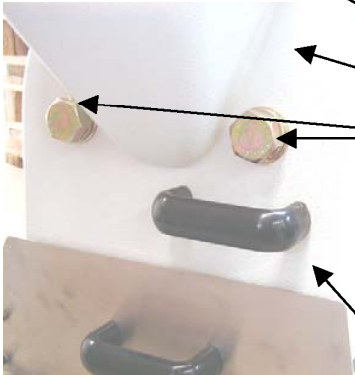

Replace the Milling Disc in the Housing:

- pull the mains plug
- remove the collecting container
- open the door
- fold up the cover
- unscrew the hexagonal screws with a SW 30 open-jawed spanner
- remove the milling disc
- clean the receptacle of the milling disc
- attach a new milling disc
- screw in the hexagonal screws again
-

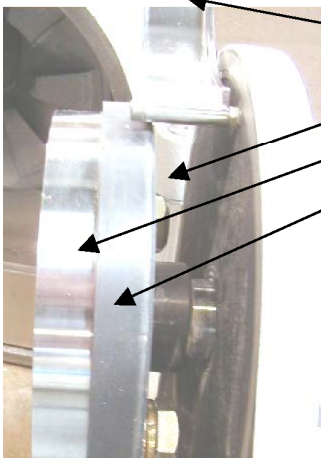

Attention !

Tighten zirconium oxide grinding disk only hand-tight – spring washer can be compressed.

Tighten steel and tungsten carbide grinding disks with a tightening torque of 430 Nm

General Information

Accessories

Cleaning

You may pull out the collecting container for cleaning purposes. Open the door to clean the milling discs, the milling space and the funnel tube.

You may now conveniently clean the entire area that has come into contact with the material to be crushed with the aid of a brush or vacuum cleaner.



Do not clean the DM 200 with flowing water.

Danger from Current Surge

Only use moistened rags.

Solvents are not permitted.

The protection type of **DM 200** is IP 55



Service

The DM 200 does not require service, although it should be relubricated on the lubricating nipple of the shifting sled.

Conventional lubricating grease may be used for this purpose.

A grease gun is also available as an accessory under the order no. 05.185.0006.

Necessary Safety Tests

The safety limit switches on the door and the collecting container should be checked on a monthly basis to see if they are in perfect working order:

- open the door:

It may not be possible to switch on the DM 200 any more with the ON/OFF switch.

- close the door:

It must be possible to start the DM 200 with the ON/OFF switch again.

- pull out the collecting container:

It may not be possible to switch on the DM 200 any more with the ON/OFF switch.

- push in the collecting container:

It must be possible to start the DM 200 with the ON/OFF switch again.



If the machine is not functioning properly, do not continue to operate it and get in contact with our customer service.

Copyright

It is only allowed to pass on or copy this documentation or utilise and pass on its content with the express permission of Retsch GmbH & Co.KG.

Non-compliance engages the obligation of compensation for damage.

Modifications

We reserve the right to make technical modifications.

Safety Regulations (Table)

for the **DM 200** from the sections

process	action	hazards
safety information	property damage and personal injury may be incurred if the safety information is not observed.	no claims to compensation for damage shall be accepted whatsoever
intended use	do not mill any explosive self-igniting or fire-promoting materials	danger to life and limb because of explosion because the DM 200 is not designed to be explosion protected.
	do not mill materials whose fine proportion at a certain percentage may lead to an explosion.	danger to life and limb because of explosion because the DM 200 is not designed to be explosion protected.
	do not modify the machine and only use the spare parts and accessories approved by Retsch.	otherwise the conformity with European directives declared by Retsch shall be rendered invalid. Furthermore, this leads to the loss of all warranty claims.
packaging	please store the packaging for the warranty period.	if there is a complaint and it is returned in insufficient packaging, this may jeopardise your warranty claims.
temperature fluctuations	protect the DM 200 from condensate water forming if there are temperature fluctuations.	electronic component parts may be damaged.
transportation	do not shock, shake or throw the DM 200 during transportation	electronic and mechanical component parts may be damaged.
scope of delivery	if the delivery is incomplete and / or if there is damage from transportation, you have to inform the transporter and Retsch GmbH without delay (within 24h).	under certain circumstances later complaints may not be recognised.
ambient temperature	should not fall below 5°C should not go above 40°C	electronic and mechanical component parts may be damaged and performance data may change unpredictably.
humidity	do not exceed a relative humidity of 80% at 31°C or 50% (descending linearly) at 40°C.	higher humidity may damage the electronic and mechanical component parts and performance data may change unpredictably.
electrical connection	the network does not agree with the values on the nameplate.	electronic and mechanical component parts may be damaged.
adjusting the gap width	do not adjust a gap width < 0.1 mm	if the milling discs come into contact with one another, they may be damaged.
feeding the material to be crushed	only feed material to be crushed when the DM 200 is operating.	if the DM 200 starts up filled with material to be crushed, this may damage the mechanical component parts.
cleaning	do not clean the DM 200 with flowing water.	danger from current surge. the protection type of the DM 200 is IP55
safety tests	get in contact with customer service if the safety equipment functions improperly..	if defects are improperly repaired, this may lead to danger to life and limb.

The Terms of the Warranty

1. We will remedy defects or furnish replacement at no charge if complaints are justified.

The purchaser shall only be entitled to repudiation of contract or diminution if according to our decision it is not possible to remedy the defect or replacement cannot be made or the period for this cannot be maintained or it was our fault that an appropriate subsequent period set by the customer of at least six weeks was not maintained.

If the remedy or supplying the replacement is finally unsuccessful, the customer may demand that the remuneration be reduced or withdraw from the contract at his choice. Further claims, especially to compensation for damage that were not incurred on the object itself, such as a loss of production, shall be excluded to the extent that we are not guilty of intent or negligence. We pass on the manufacturer's liability for external products.

2. We shall bear the immediate costs incurred by remedying the defect or supplying a replacement under the proviso that the complaint has turned out to be justified. This shall also apply with reference to the costs of shipping and the appropriate costs of removing and installing. However, the customer shall bear the appropriate costs for providing his own installers and help on site. To the extent that our customer works abroad, we shall be entitled as a deviation from this to pay the costs necessary for remedying the defect, especially transportation, travelling and material costs from the German border.

3. Our warranty is always six months. We do not assume any warranty with wear parts.

If the user fills out the warranty card attached to the goods and sends it to us within two weeks of purchasing the machine, our warranty period shall be increased to two years. We do not assume any warranty with wear parts. The warranty period shall be calculated from the time of sending the warranty card if putting the machine into service later is reported in the warranty card from the time of putting the machine into service. The machine may not be put into service any later than one half year after purchasing the machine.

The warranty refers to its use in a laboratory in 1-shift operation. The warranty period shall be reduced accordingly for multiple-shift operation or other areas of application.

4. We shall assume the warranty for our goods being free of manufacturing deficiencies. The suitability, classification and functioning of our goods shall be exclusively determined according to the performance specifications in the order confirmation, even if this deviates from the order. In this case, the customer has the possibility of pointing out any differences to the order within two weeks after receiving our order confirmation and reaching an agreement with us on this. If he does not object to the specifications in the order confirmation, it shall be deemed as accepted.

In the absence of a deviating agreement, we shall not be liable for the article of sale being suited to the purpose planned by the customer. The shall apply to the performance data that the customer expects unless we have been able to make appropriate practical laboratory experiments ahead of time and declared the corresponding performance data to be binding in our order confirmation in writing.

5. Our warranty shall also be invalid if persons other than those commissioned by us do repairs or otherwise interfere with or modify the goods supplied by us or do not use the suitable accessories to the extent that the deficiency occurring is in causal connection with that. Otherwise, the prerequisite for our warranty shall be complying with our operating instructions and instructions for use.

6. If the customer builds the goods in other systems or production plants or connects them up to them, or interlinks or processes them without our prior approval, our warranty shall be exclusively limited to the components supplied by us.

7. We may either remedy defects or exchange defective components at our choice at the place of installation of the object of sale or in our company. If the defect is remedied at the place of installation, the customer has to grant our representative temporally and spatially unobstructed access to the object of sale. Otherwise, the customer may only demand that the warranty work only be done during the business hours normal at the location. Should the warranty work be done outside of our normal business hours at the request of the purchaser, the customer has to bear the additional costs. If he wishes other special services going beyond the warranty work, these expenses shall be paid at the prices of us valid at the time.

Distributed by:



Carl Stuart Limited

ADVANCED APPLIED TECHNOLOGIES

Contact Us:

Irl Ph: 01 4523432

UK Ph: 08452 30 40 30

Web: www.carlstuart.com

Email: info@carlstuart.com