

Viscosity Standards

Brookfield Viscosity Standards provide a convenient, reliable way to verify the calibration of your Brookfield Laboratory Viscometer/Rheometer. Brookfield Viscosity Standards are

Newtonian, and they are available as either silicone or oil. Silicone fluids are less temperature sensitive than oil fluids.

Note: Brookfield recommends that all fluids be replaced annually

Silicone Viscosity Standards

These fluids are most commonly used to verify calibration of Brookfield Viscometers/Rheometers.

Accuracy: ±1% of viscosity value

Excellent temperature stability

Recommended for use with Brookfield and most other rotational viscometers


Most economical

Special viscosity values and temperature calibrations available upon request



VisCal Kit

The Brookfield VisCal Kit provides all the necessary items to verify calibration of your Viscometer/Rheometer. Includes Brookfield 600mL Beaker, 1 pint of Silicone Viscosity Standard, Dispersing Bottle for cleanup and Trapper Cleaning Agent.*

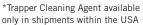



Plastic VisCal Kit

The Brookfield Plastic VisCal Kit provides all the necessary items to verify calibration of your Viscometer/ Rheometer in a glass-free environment. Includes Brookfield 600mL Plastic Beaker, 1000ml of Silicone Viscosity Standard (5-12,500 cP) in a plastic bottle and a Brookfield-designed metal lid for anchoring beaker in the temperature bath.



General Purpose Silicone Fluids				
Brookfield Part #	Nominal Viscosity cP (mPa•s)	Temp °C		
5 cps	5	25.0°C		
10 cps	10	25.0°C		
50 cps	50	25.0°C		
100 cps	100	25.0°C		
500 cps	500	25.0°C		
1000 cps	1,000	25.0°C		
5000 cps	5,000	25.0°C		
12500 cps	12,500	25.0°C		
30000 cps	30,000	25.0°C		
60000 cps	60,000	25.0°C		
100000cps	100,000	25.0°C		

High Temperature Silicone Fluids				
Brookfield Part #	Nominal Viscosity cP (mPa•s)	Temp °C	Temp °F	
HT30000	30,000	25.0°C	77°F	
	9,000	93.3°C	200°F	
	4,500	149.0°C	300°F	
HT60000	60,000	25.0°C	77°F	
	18,000	93.3°C	200°F	
	9,000	149.0°C	300°F	
HT100000	100,000	25.0°C	77°F	
	30,000	93.3°C	200°F	
	15,000	149.0°C	300°F	

Special Order Silicone Fluids

For our customers needing a nonstandard viscosity or temperature range, our silicone fluids can be modified to meet most requirements.

VISCOSITY BLENDS CALIBRATED AT 25°C (77°F)

- Minimum: 5 cP (mPa•s)

- Maximum: 60,000 cP (mPa•s)

- Blends will be within ±2% of requested value

TEMPERATURE CALIBRATIONS

- Minimum: 10°C (50°F)

Maximum: 80°C (176°F)

- Minimum temperature increment: 2°C

Distributed by: Lab Unlin Tallaght Business Park Whitestown, Dublin 24, Ireland D24 RFK3

Tel: (01) 4523432 Tel: 08452 30 40 30 Fax: 08452 30 50 30

Camberley,

GU16 7ER

United Kingdom

Fax: (01) 4523967

E-mail: info@labunlimited.com E-mail: info@labunlimited.co.uk Web: www.labunlimited.com Web: www.labunlimited.co.uk

Quatro House, Frimley Road,

Oil Viscosity Standards

These fluids are used for specific instruments using cone/plate or Krebs spindle geometry. Also, certain industries may require use of oil standards.

Accuracy: $\pm 1\%$ of viscosity value Appropriate for use at shear rates greater than 500 sec-1 Recommended for use with cone/plate Viscometers at viscosities above 5,000 cP Recommended for Brookfield CAP series and KU-2 Viscometers and RST Rheometers Brookfield oil viscosity

Note: Other oil fluids are available - call for details

standards are hydro-

carbon based, either mineral oil or polybutenes

Brookfield Viscosity Standards are accurate to ±1% of the stated viscosity and are certified by methods traceable to the United States National Institute of Standards and Technology (NIST). The selection of one or two fluids will normally provide sufficient measurement points to verify calibration of your instrument. All fluids are supplied in 1/2 liter (1 pint) containers complete with a certificate of calibration. CAP Oil Fluids are supplied in 150 mL (4 oz) containers

CAP Viscometer Oil Fluids For calibrating CAP Series cones each spindle has its own fluid								
	HIGH TORQUE CAP Low Temp 25°C High Temp 60°C			LOW TORQUE CAP Low Temp 25°C High Temp 60°C				
Cone Spindle		l Viscosity cP (mPa•s)		d Viscosity cP (mPa•s)	Brookfie Part #	eld Viscosity cP (mPa•s)		d Viscosity cP (mPa•s)
1	CAP1L	89	CAP1H	89	CAPOL	57	CAPOH	57
2	CAP2L	177	CAP2H	177	CAP1L	89	CAP1H	89
3	CAP3L	354	CAP3H	354	CAP2L	177	CAP2H	177
4	CAP4L	708	CAP4H	708	CAP3L	354	CAP3H	354
5	CAP5L	1,417	CAP5H	1,417	CAP4L	708	CAP4H	708
6	CAP6L	3,542	CAP6H	3,542	CAP5L	1,417	CAP5H	1,417
7	CAP7L	1,328	CAP7H	1,328	CAP1L	89	CAP1H	89
8	CAP8L	5,313	CAP8H	5,313	CAP3L	354	CAP3H	354
9	CAP9L	21,250	CAP9H	21,250	CAP5L	1,417	CAP5H	1,417
10	CAP10L	236	CAP10H	l 236	CAP2L	177	CAP2H	177

HOW TO SELECT A CAP FLUID

- Determine which viscometer is being used: High Torque or Low Torque.
- Determine which temperature model is being used: Low Temperature (5°C-75°C) or High Temperature (50°C-235°C)
- Determine which cone is being used.

Krebs Viscometer Oil Fluids				
Brookfield Part #	Nominal Viscosity Krebs Units	Temp °C		
KU61	61	25.0°C		
KU73	73	25.0°C		
KU87	87	25.0°C		
KU99	99	25.0°C		
KU106	106	25.0°C		

General Purpose Oil Fluids			
Brookfield Part #	Nominal Viscosity cP (mPa•s)	Temp °C	
B29	29	25.0°C	
B200	200	25.0°C	
B400	400	25.0°C	
B600	600	25.0°C	
B1060	1,060	25.0°C	
B2000	2,000	25.0°C	
B10200	10,200	25.0°C	
B21000	21,000	25.0°C	
B73000	73,000	25.0°C	
B200000	200,000	25.0°C	
B360000	360,000	25.0°C	

RST Rheometer Oil Fluids (calibrated at 25.0°C)				
Cone Spindle	Brookfield Part #	Nominal Viscosity cP (mPa•s)		
RCT-25-1	B41000	41,000		
RCT-25-2	B73000	73,000		
RCT-50-1	B10200	10,200		
RCT-50-2	B21000	21,000		
RCT-75-1	B4900	4,900		
RCT-75-2	B10200	10,200		

RST Rheometer Oil Fluids (calibrated at 25.0°C)			
Coaxial Spindle	Brookfield Part #	Nominal Viscosity cP (mPa•s)	
CCT-DG	B200	200	
CCT-40	B2000	2,000	
CCT-25	B10200	10,200	
CCT-14	B73000	73,000	
CCT-8	B360000	360,000	