PRECISION



SADE SP140

Weight Sorter/Checkweigher for tablets and capsules

When precision and accuracy count

The SADE SP140 is ideal for both Research & Development and Production departments for the weight sorting of batches of tablets and capsules. Each item is weighed to the high accuracy of +/- 1mg and sorted into accept and reject categories.

It is a floor standing unit with one weighing channel automatically fed by a 60 litre bulk hopper. Once set up the SADE SP140 can be left unattended for long periods of operation freeing operators for other tasks. Dismantling for cleaning is rapid and it's easy to set up for the next product

R&D — Ensures weight uniformity

Lab Unlimited

Distributed by

100% weight checking and statistical analysis of clinical trial or development batches. SADE SP units are also ideal for checking over encapsulation.

PRODUCTION — Recover valuable product

Enables you to reclaim your weight acceptable product from rejected batches – reducing wastage and cost.

Savings can be made recovering valuable product made during press run-in and run-out stages.

STILART

Tallaght Business Park

Fax: (01) 4523967 E-mail: info@labunlimited.com

Web: www.labunlimited.com

Whitestov Dublin 24

Tel: (01) 4523432

Ireland

KEY FEATURES

- Sorts both tablets and capsules
- Handles different shapes and sizes
- No change parts required
- High accuracy, +/- 1mg
- Quick and easy clean down
- 24 hour operation
- Minimal operator involvement



SADE SP140

Weight Sorter/Checkweigher for tablets and capsules

KEY OPERATIONAL BENEFITS

Our SADE SP units are already delivering precision weight sorting of the world's leading pharmaceutical companies around the globe.


Maximise speed and accuracy, minimise cost

- Consistent and precision weight sorting with an accuracy of +/- 1mg
- Versatile unit sorts tablets and capsules of different shapes and sizes
- No need for costly format parts and time consuming change over's
- Fast effective clean down boosts operational efficiency
- Quick and easy batch set up via improved user interface
- Minimal operator involvement required saving time and funds
- Proven design and quality materials for ease of cleaning

Informed decision making and action taking

- New touch screen provides real time display of statistical data
- Prints detailed reports directly via USB local or network printer
- Data can be transferred to your PC via SP Connect
 - enhanced reports automatically saved as secure pdf files
 - individual weights can be transferred into Excel



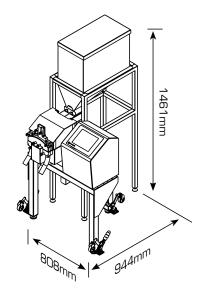

TECHNICAL DETAILS

Floor standing
1
Up to 4500 items per/hr Varies per product and weighing environment
Up to 105,000 per 24 hour period
1 x 60 Litres
0 - 2g
Tablets: Up to 24mm length, 20mm wide Capsules: 0,1,2,3,4,5,& 00
0.1 mg
+/- 1mg
Automatic
316 electro-polished stainless steel and PEEK polyetheretheketone FDA 21 CFR 177.2415
944D x 808W x 1461H
115/230 VAC 50/60 Hz
Clean dry air 3-5 bar or 43-73PSI



PRECISION

© CI Precision 2011 — SP140 Eur 0711-1.0