

# Data sheet

## dGTP

(2'-deoxy-guanosine-5'-triphosphate)

**Cat. No N0095**

40  $\mu\text{mol}$  (0.4 mL of 100 mM solution)

### Description

Ultrapure dGTP supplied as clear aqueous solution (pH 8.5).


**Formula:**  $\text{C}_{10}\text{H}_{16}\text{N}_5\text{O}_{13}\text{P}_3$

**Molecular Weight:** 507.18 g/mol

$\lambda_{\text{max}}$  pH 7.0 = 252 nm

$\epsilon$  at  $\lambda_{\text{max}}$ , pH 7.0 = 13.7  $\text{mmol}^{-1}\text{cm}^{-1}$

**Purity:** > 99% confirmed by HPLC

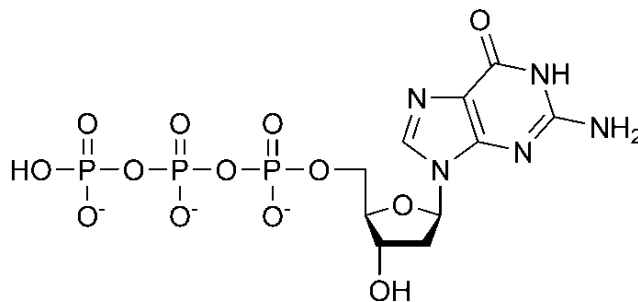

### Applications

For use in all molecular biology applications, including PCR, real-time PCR, high fidelity and long PCR, cDNA synthesis, RT-PCR, and DNA sequencing.

The synthesis of a complementary DNA/cDNA strand by enzymatic polymerization of deoxynucleotide triphosphates (dNTPs) is the basic principle of DNA amplification, sequencing and cDNA synthesis techniques.

**Concentration:** 100 mM

### Quality Control Specifications:

18 kb long range PCR (template dilution series): suitable; contamination with bacterial and human DNA: not detectable; activity of DNase, Protease or Phosphatase: not detectable.

**Storage:** Store at  $-20^{\circ}\text{C}$ , short term (up to one week) exposure to ambient temperature is possible.

### PRODUCT USE LIMITATION

This product is developed, designed and sold exclusively for research purposes and in vitro use only. The product was not tested for use in diagnostics or for drug development, nor is it suitable for administration to humans or animals. Please refer to [www.canvaxbiotech.com](http://www.canvaxbiotech.com) for Material Safety Data Sheet of the product.

Distributed by:

**Lab Unlimited**  
CARL STUART GROUP

Tallaght Business Park  
Whitestown, Dublin 24,  
Ireland  
D24 RFK3

Tel: (01) 4523432  
Fax: (01) 4523967  
E-mail: [info@labunlimited.com](mailto:info@labunlimited.com)  
Web: [www.labunlimited.com](http://www.labunlimited.com)

Quatro House, Frimley Road,  
Camberley,  
United Kingdom  
GU16 7ER

Tel: 08452 30 40 30  
Fax: 08452 30 50 30  
E-mail: [info@labunlimited.co.uk](mailto:info@labunlimited.co.uk)  
Web: [www.labunlimited.co.uk](http://www.labunlimited.co.uk)

 **canvax**

## Notes

Distributed by:

**Lab Unlimited**  
CARL STUART GROUP

Tallaght Business Park  
Whitestown, Dublin 24,  
Ireland  
D24 RFK3

Tel: (01) 4523432  
Fax: (01) 4523967  
E-mail: [info@labunlimited.com](mailto:info@labunlimited.com)  
Web: [www.labunlimited.com](http://www.labunlimited.com)

Quatro House, Frimley Road,  
Camberley,  
United Kingdom  
GU16 7ER

Tel: 08452 30 40 30  
Fax: 08452 30 50 30  
E-mail: [info@labunlimited.co.uk](mailto:info@labunlimited.co.uk)  
Web: [www.labunlimited.co.uk](http://www.labunlimited.co.uk)