

NEW DV1™ Viscometer

The only viscometer in its class to offer continuous sensing and data display at such a low price!

User Configuration Display

- User choice of most important parameter is displayed in larger font size
- Choice of static or scrolling display mode

Displayed Info:

- Viscosity (cP, P, mPa•s, Pa•s)
- % Torque
- Speed/Spindle
- Temperature (°C or °F)
if RTD Temperature Probe is purchased

Choice of Multiple Languages

English, French, German, Portuguese, Russian, Spanish

USB PC interface for use with optional Wingather SQ Software

18 speeds

provide great range capability

Optional RTD Temperature Probe DVP-94Y

Direct access

to time measurement function (time to torque, time to temperature, time to stop)

Accuracy:

±1.0% of range

Repeatability: ±0.2%

Temperature off-set

capability to ±5°C

Automatic Range Calculation:

- Full Scale Range (FSR) at 100%
- Maximum viscosity measured with Spindle/Speed combination

Simplified User interface

for more direct access to features

Printing to Dymo® Printer Capability

Stylish Model G Base



What's Included?

Instrument
6 spindles (RV/HA/HB)
or 4 spindles (LV)
Spindle Guard Leg*
Lab Stand (Model G)
Carrying Case

*Not applicable to HA or HB torque models

Optional Accessories

Wingather SQ Software ►
RTD Temperature Probe
Ball Bearing Suspension (RV/HA/HB)
EZ-Lock Spindle Coupling System
Viscosity Standards
Protective Keypad Covers
Dymo Printer
RV/HA/HB-1 Spindle
Quick Action Lab Stand
Temperature Bath
Small Sample Adapter
UL Adapter
Thermosel
Helipath Stand with T-bar Spindles
Spiral Adapter
DIN Adapter
Quick Connect/Extension Links
Vane Spindles

MODEL	VISCOSITY RANGE cP(mPa•s)		SPEEDS	
	Min.	Max.	RPM	Number of Increments
DV1MLV	1*	2M	.3-100	18
DV1MRV	100	13M	.3-100	18
DV1MHA	200	26M	.3-100	18
DV1MHB	800	104M	.3-100	18

* Minimum ranges can be extended to as low as 1 cP with the use of Brookfield Accessories

** Standard torque range values

M = 1 million cP = Centipoise mPa•s = MilliPascal•seconds

Wingather SQ Software Optional (see p14 for more details)

DATA COLLECTION SOFTWARE TO COLLECT, ANALYZE AND RECORD TEST DATA

Wingather software provides an easy way to gather data and plot graphs while creating permanent test records.

Important features and benefits enhance operator versatility in performing viscosity tests:

- Multiple test modes to enhance data collection
- Follow up events including analysis through math models which calculate yield stress and plastic index
- Automatic sample numbering
- Data graphing of up to 20 data sets concurrently
- Data export to spread sheet format (Excel®)




The DV1 Viscometer communicates to the PC through USB A port. The interface cable is supplied with Wingather SQ Software. Successful communication is indicated by a green light beside the port designation.




DV1CP

The DV1 is available in a Wells/Brookfield Cone & Plate Version
Must be ordered when instrument is first purchased.

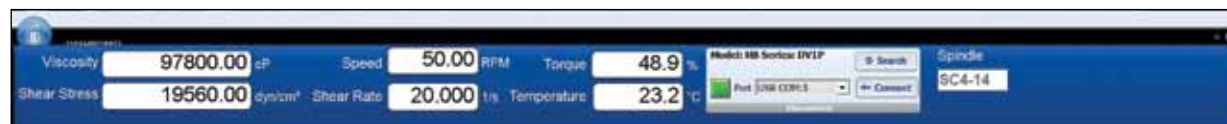

Dashboard shows current test information.

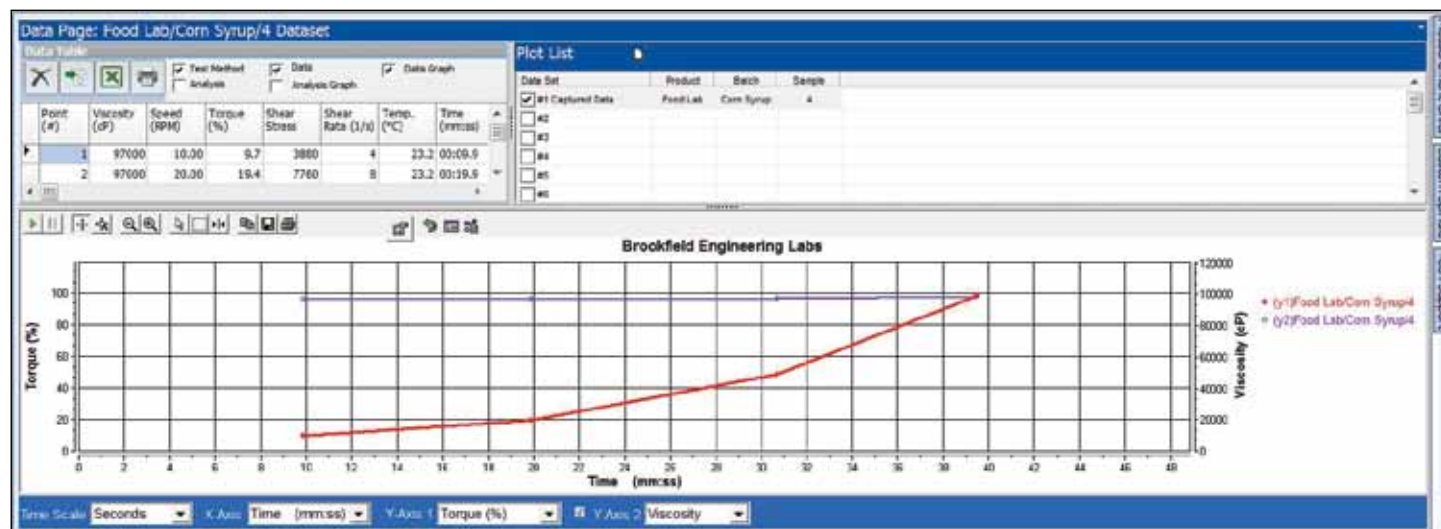

Data table and graph shows test record at a glance.

Excel® is a registered trademark of Microsoft Corporation in the United States and/or other countries.



Distributed by:

Lab Unlimited
CARL STUART GROUP

Tallaght Business Park
Whitestown, Dublin 24,
Ireland
D24 RFK3

Tel: (01) 4523432
Fax: (01) 4523967
E-mail: info@labunlimited.com
Web: www.labunlimited.com

Quatro House, Frimley Road,
Camberley,
United Kingdom
GU16 7ER

Tel: 08452 30 40 30
Fax: 08452 30 50 30
E-mail: info@labunlimited.co.uk
Web: www.labunlimited.co.uk