


Lovibond® Water Testing Tintometer® Group



TC 140 G

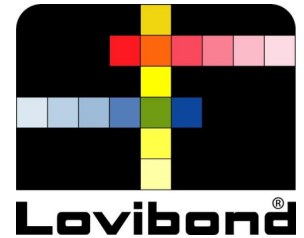

- Temperature range 2 °C to 40 °C, adjustable in steps of 0.1 °C
- Low power consumption
- Illuminated LED display of preset and current temperatures
- Ideal for BOD determination at 20 °C

Part Number: 2438210

Versatile in use

Our thermostatically controlled cabinets are the means of choice both for temperature-controlled storage of samples and for BOD determination in waste water analysis. The integral temperature control unit meets EMC directive IEC 61326 ("Electrical devices for measurement, monitoring and for use in laboratories").

Continuous temperature control

The thermostat cabinets allow continuous temperature control. Over a range of 2-40 °C the temperature can be set in steps of 0.1 °C.

Internal power connections

The use of external devices is no problem thanks to several sockets installed in the interior.

Energy efficiency

The highly efficient components and reinforced insulation provide low consumption values.

Variable

You can choose between 4 models with steel door (135, 175, 255 and 445 l net capacity) and 2 models with glass door (140 and 235 l net capacity).

Industry

Chemical Industry | Food and Beverage Industry | Industries Others | Marine Industry | Municipalities | Oil and Gas | Pharmaceutical Industry | Power and Energy

TC 140 G

Perfect for a number of different applications in industry and research. The temperature of the thermostatically controlled cabinets can be regulated individually in the smallest of steps - and an LED display keeps the status quo under control at all times. There are six models available with different useful contents and glass or steel doors.

Technical Data

Display	Backlit LED display
Power Supply	220 - 240 V / 50 Hz
Net Capacity	140 L
Consumption	approx. 1,61 kWh/ 24h
Climate Class	SN, + 10 °C to + 32 °C
Coolant	R600a
Sockets	4
Temperature Range	+ 2 °C to + 40 °C, steps of 0,1 °C
Temperature Tolerance	For BOD (T = 20 °C ±1 °C according to EN 1899)
Interior Light	LED light strip
Storage Space	3 metal grids + 1 bottom grid
Shelf Loading Capacity	45 kg
Suitable to build under	Yes
Further Equipment	Standard Door

Compliance	CE
Dimensions	600 x 850 x 600 mm
Dimensions Internal	513 x 702 x 441 mm
Weight	48 kg

Accessories

Title	Part Number
Temperature control unit for TC 135 S, TC 140 G	2428201-2
Wooden box for TC 135 S, TC 140 G (IPPC Standard - ISPM15)	425768
Metal grid for TC 135 S, TC 140 G	425791
Spare board for temperature control unit	427542

Tintometer GmbH
 Lovibond® Water Testing
 Schleefstraße 8-12
 44287 Dortmund
 Tel.: +49 (0)231/94510-0
 sales@lovibond.com
 www.lovibond.com
 Germany

The Tintometer Limited
 Lovibond House
 Sun Rise Way
 Amesbury, SP4 7GR
 Tel.: +44 (0)1980 664800
 Fax: +44 (0)1980 625412
 sales@lovibond.uk
 www.lovibond.com
 UK

Tintometer China
 Room 1001, China Life Tower
 16 Chaoyangmenwai Avenue,
 Beijing, 100020
 Tel.: +86 10 85251111 App. 330
 Fax: +86 10 85251001
 chinaoffice@tintometer.com
 www.lovibond.com
 China

Tintometer South East Asia
 Unit B-3-12, BBT One Boulevard,
 Lebuhr Nilam 2, Bandar Bukit Tinggi,
 Klang, 41200, Selangor D.E
 Tel.: +60 (0)3 3325 2285/6
 Fax: +60 (0)3 3325 2287
 lovibond.asia@tintometer.com
 www.lovibond.com
 Malaysia

Tintometer Brazil
 Caixa Postal: 271
 CEP: 13201-970
 Jundiaí – SP
 Tel.: +55 (11) 3230-6410
 sales@tintometer.com.br
 www.lovibond.com.br
 Brazil

Tintometer Inc.
 6456 Parkland Drive
 Sarasota, FL 34243
 Tel: 941.756.6410
 Fax: 941.727.9654
 sales@lovibond.us
 www.lovibond.com
 USA

Tintometer India Pvt. Ltd.
 Door No: 7-2-C-14, 2nd, 3rd & 4th Floor
 Sanathnagar Industrial Estate,
 Hyderabad, 500018
 Telangana
 Tel: +91 (0) 40 23883300
 Toll Free: 1 800 599 3891/ 3892
 indiaoffice@lovibond.in
 www.lovibondwater.in
 India

Tintometer Spain
 Postbox: 24047
 08080 Barcelona
 Tel.: +34 661 606 770
 sales@tintometer.es
 www.lovibond.com
 Spain

Technical changes without notice
 Printed in Germany

Lovibond® and Tintometer® are Trademarks of the Tintometer Group of Companies